
Preemie Size 1
& Size 2

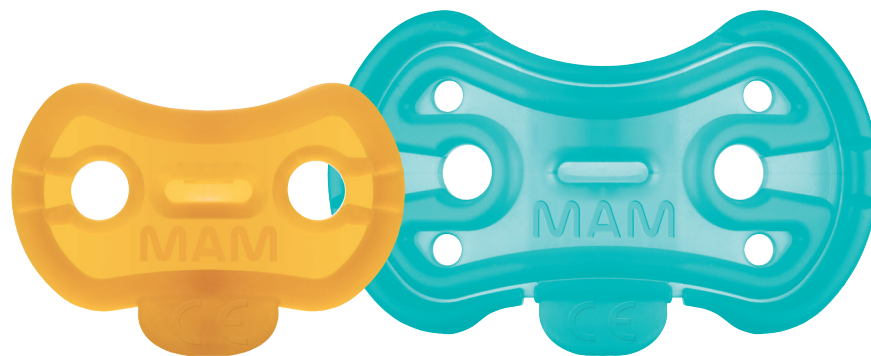Medical Device Class 1

Size 1: < 1,000 g
Size 2: 900 – 1,500 g

Select your country

Austria / Germany  

Belgium / Netherlands  

Cyprus 

Finland 

Greece 

Sweden 

Switzerland 

UK 

US 

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
5. For the child's safety – WARNING! ⚠

WARNING – This Pacifier Is For A Premature Infant. Use Of This Product Beyond The Recommended Age For The Baby Could Pose A Choking Risk.

Do Not Use This Product When the Infant Reaches the Developmental Age of 27 Weeks (<1,000 g) in case of Preemie 1 or 30 weeks (<1,500 g) in case of Preemie 2. Please Consult With Your Doctor.

Inspect carefully before each use. Pull the soother in all directions. Throw away at the first signs of damage or weakness. Do not leave the soother in direct sunlight or near a source of heat (except during cleaning) or leave in cleansing solution for longer than recommended, as this may weaken the teat. To avoid danger of suffocation, keep the plastic bag away from babies and children. Do not use the bag in cribs, beds, carriages or playpens. The bag is not a toy. Do not tie soother around child's neck as it presents a strangulation danger. The product does not meet the safety standards which are applied to normal soothers, namely EN 1400 and 16 CFR 1511. It poses a choking hazard. Therefore, it must only be used under medical supervision. Parents must never be allowed to take the product home.

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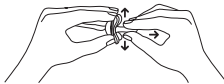5°C/ 41°F – 30°C/ 86°F The soother has to be stored at room temperature between +5°C and +30°C / +41°F and +86°F.

Protect from direct sunlight.

Expiry date is on the packaging. Replace the soother at least every month.
Dispose of this product in general hospital waste.

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1. The benefits of soothers


MAM Preemie soothers are intended to be used for premature babies up to a maximum body weight of 1,500 g under medical supervision to compensate the typically present suck and swallow deficiency. Soothers help to satisfy the natural sucking urge of babies. Babies can benefit from the pain-relieving effect of soothers during medical interventions. Furthermore, soothers help with the switch from gavage feeding to breastfeeding or bottle feeding as it trains the oral musculature. The soother is one of the many preventive strategies against sudden infant death syndrome (SIDS).

2. When to use this product

The MAM Preemies are intended for premature babies up to 1,500g under medical supervision. MAM Preemie size 1 for babies with a weight below 1,000g MAM Preemie size 2 for babies with a weight between 900g and 1,500g The soother is only intended for a premature baby in the respective weight range (see print on bag) and may only be used during hospital stay.

3. How to clean & use this product

Germ-reduced through gamma irradiation and ready for initial use if package is not opened or damaged. If packaging is damaged, clean and autoclave before use. Allow to cool. Clean after each use. Wash with mild soapy water and then rinse with clear water. Do not use aggressive cleaning agents. Autoclavable up to to 134°C/273°F for a maximum of 18 minutes (max. 20 times). Allow to stand for at least 5 minutes before giving it to baby. After cleaning or autoclaving expel remaining water by squeezing the teat with sterilised teat tongs. Replace the soother at least every month. Only use for one baby. Sizes and weights are general guidelines. Individual adaption and usage with oral & nasal tubing: The product can be tailored by cutting with scissors or a scalpel (see fig. 1 and fig. 2). Dispose of all small parts immediately after cutting to prevent choking hazards.

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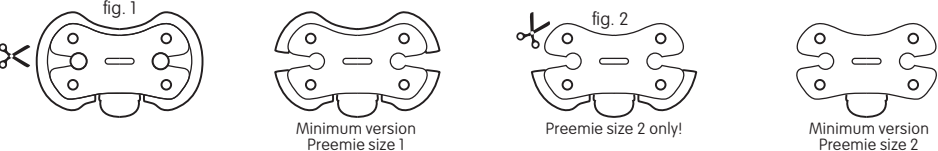MAM Premies: Perfected in cooperation with clinical experts.

Specially designed for the needs of premature babies.

7 years of intensive research and development have gone into MAM Premies 1 and 2. The soothers for premature babies have been perfected together with experts such as dentists, paediatricians, paediatric nurses, speech therapists, engineers and designers. The result is a **one-piece construction made of 100 % silicone**, which meets the special requirements for everyday use in hospitals. The MAM Premies are available in 2 sizes. The MAM Premies 1 & 2 have been developed for premature babies weighing less than 1,000 g to 1,500 g. The MAM Comfort completes the MAM Premie range and is designed for babies of 1,400g and above. The MAM Comfort can be taken home by parents.

Clinical experts recognise the positive features:

MAM Premies can be easily held in the mouth due to their low weight, and so it is not necessary to repeatedly offer them to the baby.

The **small shield design** leaves more space for breathing masks and other medical devices (Fig. 1). The **low weight and small size** mean that the lip and tongue muscles are not overworked (Fig. 2).

Fig. 1

Fig. 2

DEVELOPED WITH MEDICAL EXPERTS


Teamwork with medical experts for maximum safety

Right from the beginning, MAM products are developed in cooperation with medical experts and research institutes such as ICMRS, paediatricians, child development experts, dentists and midwives. Only after approval by medical experts is a MAM innovation ready for baby life.

MAM Premies & MAM Comfort™

Because every baby is special: MAM has developed a soother line especially for the needs of premature and low-birth-weight babies.

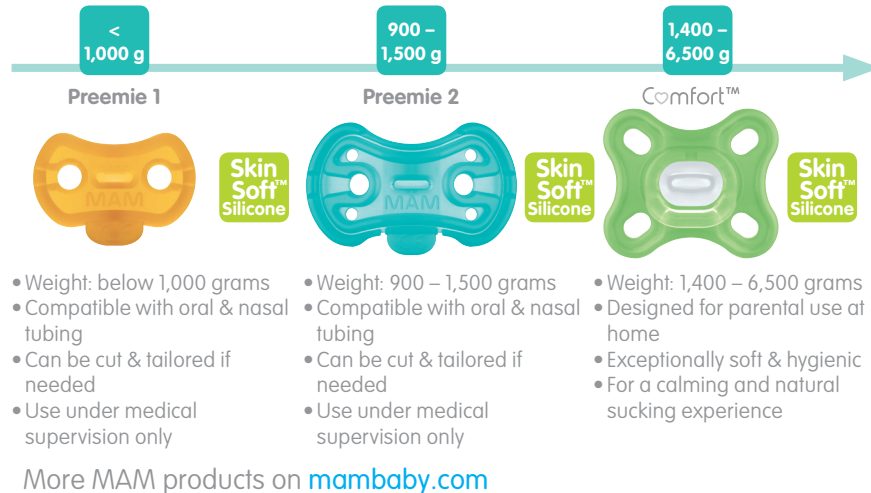For more information please contact:
MAM (UK) LTD, PO BOX 271 TEDDINGTON, TW11 8FT
Tel: 0208 943 8880
Made in Austria
Please keep instructions for future use.

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Any serious incident involving the device must be reported to the manufacturer and the competent authority of the Member State where the user and/or patient is established.

GGA1F5AC003.AAB


mambaby.com
Made in Europe

Preemie Size 1 & Size 2

- Soother for premature & low-birth-weight babies
- Ready to use – for distribution & use in hospitals only
- Lightweight one-piece construction
- Latex-free

*BPA/BPS free: All MAM products are made from materials free of BPA and BPS.

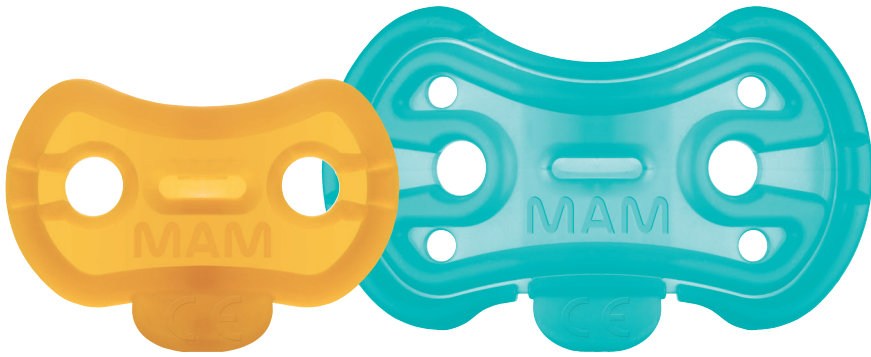Medical Device
Class I

Instructions for use – MAM Preemies size 1 & size 2

 Please read these instructions carefully before first use! These instructions only apply to MAM Preemies size 1 and size 2. Separate instructions are available for the MAM Comfort.

1. The benefits of soothers


MAM Premie soothers are intended to be used for premature babies up to a maximum body weight of 1,500 g under medical supervision to compensate the typically present suck and swallow deficiency. Soothers help to satisfy the natural sucking urge of babies. Babies can benefit from the pain-relieving effect of soothers during medical interventions. Furthermore, soothers help with the switch from gavage feeding to breastfeeding or bottle feeding as it trains the oral musculature. The soother is one of the many preventive strategies against sudden infant death syndrome (SIDS).

2. When to use this product

The MAM Preemies are intended for premature babies up to 1,500g under medical supervision.
MAM Premie size 1 for babies with a weight below 1,000g
MAM Premie size 2 for babies with a weight between 900g and 1,500g
The soother is only intended for a premature baby in the respective weight range (see print on bag) and may only be used during hospital stay.

3. How to clean & use this product

Germ-reduced through gamma irradiation and ready for initial use if package is not opened or damaged. If packaging is damaged, clean and autoclave before use. Allow to cool. Clean after each use. Wash with mild soapy water and then rinse with clear water. Do not use aggressive cleaning agents. Autoclavable from 121°C/250°F up to 134°C/273°F, for 15 (min.) to 18 (max.) minutes (max. 20 times). Allow to stand for at least 5 minutes before giving it to baby. After cleaning or autoclaving expel remaining water by squeezing the teat with sterilised teat tongs. Replace the soother at least every month. Only use for one baby. Sizes and weights are general guidelines. Individual adaption and usage with oral & nasal tubing: The product can be tailored by cutting with sterilised scissors or a sterilised scalpel (see fig. 1 and fig. 2). Dispose of all small parts immediately after cutting to prevent choking hazards.

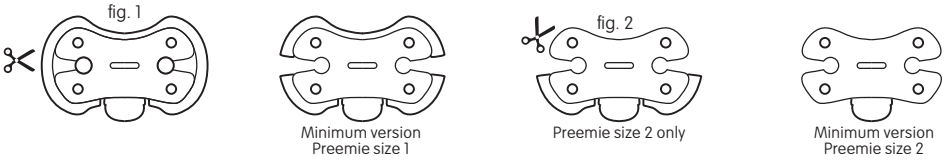4. When this product may not be used


The MAM Preemies size 1 & 2 may only be used during hospital stay under medical supervision. Therefore, parents must NEVER be allowed to take the product home.

MAM Premie size 1 must not be used for babies above 1,000g.
MAM Premie size 2 must not be used for babies below 900g or above 1,500g.

The size of the soother is embossed on the teat.

5. For the child’s safety – WARNING!

WARNING – This Soother Is For A Premature Infant. Use Of This Product Beyond The Recommended Age For The Baby Could Pose A Choking Risk.
Do Not Use This Product When the Infant Reaches the Developmental Age of 27 Weeks (<1,000 g) in case of Premie 1 or 30 weeks (<1,500 g) in case of Premie 2.
Please Consult With Your Doctor.

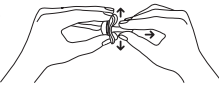Inspect carefully before each use. Pull the soother in all directions. Throw away at the first signs of damage or weakness. Do not leave the soother in direct sunlight or near a source of heat (except during cleaning) or leave in cleansing solution for longer than recommended, as this may weaken the teat. To avoid danger of suffocation, keep the plastic bag away from babies and children. Do not use the bag in cribs, beds, carriages or playpens. The bag is not a toy. Do not tie soother around child’s neck as it presents a strangulation danger. The product does not meet the safety standards which are applied to normal soothers, namely EN 1400 and 16 CFR 1511. It poses a choking hazard. Therefore, it must only be used under medical supervision. Parents must never be allowed to take the product home.

6. Composition of this product

The MAM Preemies are made of 100% odourless silicone.

7. How to store and dispose of this product

 5°C/ 41°F – 30°C/ 86°F The soother has to be stored at room temperature between +5°C and +30°C / +41°F and +86°F.

	Protect from direct sunlight.
	Expiry date is on the packaging. Replace the soother at least every month. Dispose of this product in general hospital waste.
	Lot number is on the packaging.
	Date of manufacture is on the packaging
	MAM Baby AG, Sihleggstrasse 15, 8532 Wollerau, Switzerland
	MAM Health & Innovation GmbH, Innovationsplatz 1, 7051 Großhöflein, Austria
	single patient multiple use
Date of information: 12/2020 – V.1.1	

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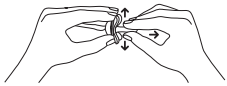
More MAM products on mambaby.com

Patent number: US 7789894

For more information please contact:

MAM USA Corporation
2700 Westchester Avenue, Suite 315
Purchase, NY 10577
914-468-4048
NICUcustomerservice@mambaby.com
Distributed by MAM BABY AG, Sihleggstrasse 15, 8832 Wollerau, Switzerland
Made in Austria
Please keep instructions for future use.
Conserve las instrucciones para uso futuro.

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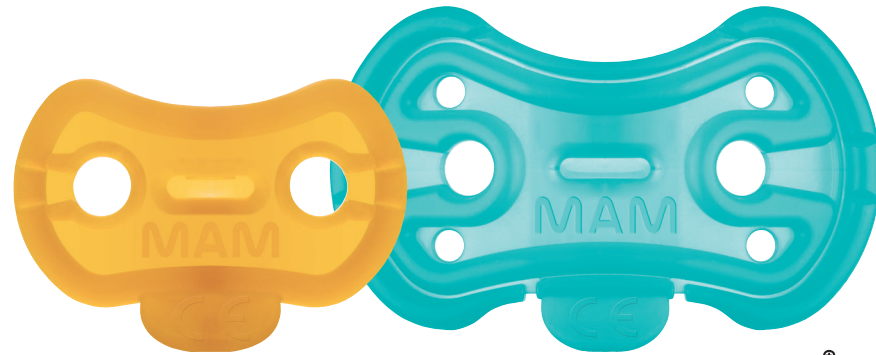
MAM Preemies size 1 & size 2

Instructions for use:

Please read these instructions carefully before first use! These instructions only apply to MAM Preemies size 1 and size 2. Separate instructions are available for the MAM Comfort.

The benefits of pacifiers

MAM Preemie pacifiers are intended to be used for premature babies up to a maximum body weight of 1,500 g under medical supervision to compensate for the typically present suck and swallow deficiency. Pacifiers help to satisfy the natural sucking urge of babies. Babies can benefit from the pain-relieving effect of pacifiers during medical interventions. Furthermore, pacifiers help with the switch from gavage feeding to breastfeeding or bottle feeding as it trains the oral musculature. The pacifier is one of the many preventive strategies against sudden infant death syndrome (SIDS).

When to use this product

The MAM Preemies are intended for premature babies up to 1,500 g under medical supervision.


MAM Preemie size 1- for babies with a weight below 1,000 g

MAM Preemie size 2- for babies with a weight between 900 g and 1,500 g

The pacifier is only intended for a premature baby in the respective weight range (see print on bag) and may only be used during hospital stay.

How to clean & use this product

MAM preemies are germ-reduced through gamma irradiation and ready for initial use if package is not opened or damaged. If packaging is damaged, clean and autoclave before use. Allow to cool. Clean after each use. Wash with mild soapy water and then rinse with clear water. Do not use aggressive cleaning agents. Autoclavable up to 134°C/273°F for a maximum of 18 minutes (max. 20 times). Allow to stand for at least 5 minutes before giving to baby. After cleaning or autoclaving expel remaining water by squeezing the nipple with sterilized nipple tongs. Replace the pacifier at least every month. Only use for one baby. Sizes and weights are general guidelines. Individual adaption and usage with oral & nasal tubing: the product can be tailored by cutting with scissors or a scalpel (see fig. 1 and fig. 2). Dispose of all small parts immediately after cutting to prevent choking hazards.

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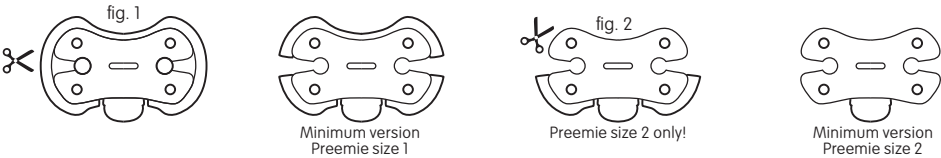
The size of the pacifier is embossed on the nipple.

⚠ WARNING! For your child's safety and health.

WARNING: This pacifier is for a premature infant. Use of this product beyond the recommended age for the baby could pose a choking risk.
Do not use this product when the infant reaches the developmental age of 27 Weeks (<1,000 g) in case of Preemie 1 or 30 weeks (<1,500 g) in case of Preemie 2. Please consult with your doctor.

Inspect carefully before each use. Pull the pacifier in all directions. Throw away at the first signs of damage or weakness. Do not leave the pacifier in direct sunlight or near a source of heat (except during cleaning) or leave in cleansing solution for longer than recommended, as this may weaken the nipple. To avoid danger of suffocation, keep the plastic bag away from babies and children. Do not use the bag in cribs, beds, carriages or playpens. The bag is not a toy. Do not tie pacifier around child's neck as it presents a strangulation danger. The product does not meet the safety standards which are applied to normal pacifiers, namely EN 1400 and 16 CFR 1511. It poses a choking hazard. Therefore, it must only be used under medical supervision. Parents must never be allowed to take the product home.

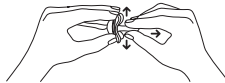Composition of this product

The MAM Preemies are made of 100% odorless silicone.

How to store and dispose of this product

5°C / 41°F 30°C / 86°F The pacifier has to be stored at room temperature between + 5°C and + 30°C / + 41°F and + 86°F.

- Protect from direct sunlight.
- Best used by date is on the packaging. Replace the pacifier at least every month. Dispose of this product in general hospital waste.
- Lot number is on the packaging.
- Date of information: 10/2018 – V.1.0

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